

CPV/CCV/Giardia Ag

General Information

Purpose : Detection of Canine Parvo, Corona and Giardia antigen

Principle : Immunochromatographic assay

Specimen : Canine feces

Materials provided :

- Anigen Rapid CPV/CCV/Giardia test device
- Assay diluents tube for CPV/CCV Ag (P)
- Assay diluents tube for Giardia Ag (G)
- Disposable swab
- Disposable dropper



Special Features :

- Simultaneous detection of CPV, CCV and Giardia
- One-step testing procedure : Quick & Accurate
- High sensitivity and specificity




Test Procedure


[Proper sample amount collected by a swab]

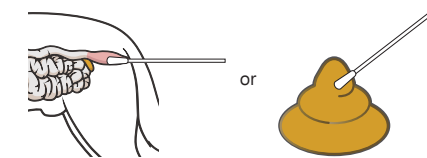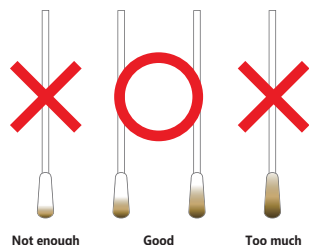


1 Collect a sample from the canine feces by using the swab

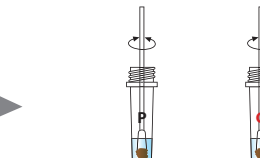


2 Insert the swab into the sample collection tube and stir the swab at least 10 times

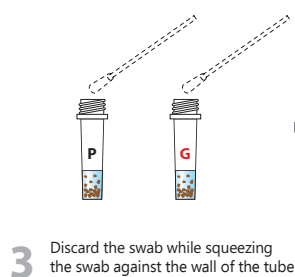

3 Discard the swab while squeezing the swab against the wall of the tube

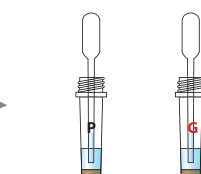


4 Wait until large particles to sink and take the supernatant by using the disposable dropper provided

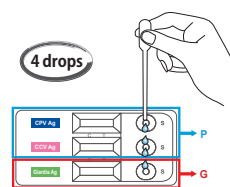

5 Add 4 drops into the sample hole on the test device

Interpretation

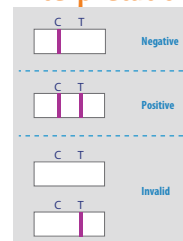

Ordering Information

Cat. No.	Description	Type	Pack size
RC1112DC	CPV/CCV/Giardia Ag	Quattro Device	5 Tests/Kit